


AQUA-CAP™ MICRODISPENSERS FROM DRUMMOND

Drummond Aqua-Cap™ Microdispensers Plug In The Accuracy

Single Use, Disposable Glass Capillary Pipets Have Hydrophobic Plug Preset to Measure Precise Volume

The Drummond Aqua-Cap preset disposable microdispensers take the guesswork out of micro liquid transfer. It is a unique, single use dispensing capillary pipet with a factory installed hydrophobic plug which determines volume. The plug controls the liquid measure and serves as a piston which is pressed downward by a reusable plunger rod. Once used, the Aqua-Cap cannot be used again, which eliminates any potential cross contamination. Available in sizes from 10 to 100 µl in packages of 100—plain or heparin coated.

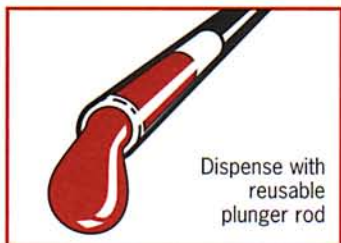

Catalog Number	Description
7-000-0100	Aqua-Cap, 10 µl, with 5 plungers
7-000-0200	Aqua-Cap, 20 µl, with 5 plungers
7-000-0250	Aqua-Cap, 25 µl, with 5 plungers
7-000-0400	Aqua-Cap, 40 µl, with 5 plungers
7-000-0447	Aqua-Cap, 44.7 µl, with 5 plungers
7-000-0500	Aqua-Cap, 50 µl, with 5 plungers
7-000-1000	Aqua-Cap, 100 µl, with 5 plungers

Mylar® clad available on special OEM orders
Custom volumes also available

Available through your local distributor



500 Parkway, Box 700
Broomall, PA 19008
Call Toll Free: 1-800/523-7480
In PA: 1-610/353-0200
FAX: 1-610/353-6204

Aqua-Cap™ is a trademark of Drummond Scientific Company.
Mylar® is a registered trademark of DuPont.
©Copyright, Drummond Scientific Company, 1996 Printed in USA Rev. 8/96